

Delphion

RESEARCH

INTEGRATED IAM

SERVICES

INSIDE DELPHION

Log Out | Work Files | Saved Searches

The Delphion Integrated View

Buy Now: [More choices...](#)Tools: [Add to Work File](#) | [Create new Work File](#)View: [INPADOC](#) | Jump to: [Top](#) | [Go to: Derwent...](#)

E

>Title: **JP3146679A2: LIVING BODY-RESTORING MEMBER MADE OF TITANIUM OR TITANIUM ALLOY AND ITS SURFACE TREATMENT**

Country: **JP Japan**
Kind: **A**

Inventor: **KAWAHARA HARUYUKI;
NOMURA YUTAKA;
TSUKAMOTO SEIICHI;**

Assignee: **KAWAHARA HARUYUKI
TOHO TITANIUM CO LTD**
[News, Profiles, Stocks and More about this company](#)



Published / Filed: **June 21, 1991 / Oct. 30, 1989**

Application Number: **JP1989000282570**

IPC Code: **C23F 1/26; A61C 8/00; A61F 2/28; A61L 27/00; B24B 1/00;**

Priority Number: **Oct. 30, 1989 JP1989000282570**


Abstract: **PURPOSE:** To make irregular fine recesses in the surface of at least the embedded part of a living body-restoring member made of Ti or Ti alloy and to enhance the adhesive strength of the member to tissue by immersing the member in an aq. HF soln. and treating it with a mixed aq. soln./contg. HF and H2O2.

CONSTITUTION: The surface of at least the embedded part of a living body- restoring member for dentistry or surgery made of Ti or Ti alloy is pretreated by immersion in a 1-6 wt. % aq. HF soln. for 30sec-3min. The pretreated surface is then treated by immersion in an aq. soln. prep'd. by mixing a 1-6% aq. HF soln. with a 1-10% aq. H2O2 soln. for 10-60sec. By this treatment, irregular fine recesses of 1-10µm average diameter and 0.5-5µm average depth are made in the surface of at least the embedded part and a certain anchoring effect can be produced between the surface of the restoring member and bone tissue, etc.

COPYRIGHT: (C)1991,JPO&Japio

Family: [Show 2 known family members](#)

Other Abstract Info: **DERABS C91-227175 DERC91-227175**

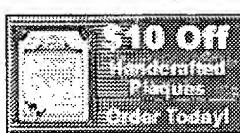

[Nominate](#)

[this for the Gallery...](#)